

# SureSeq



A Sysmex Group Company

## Application Note

# Improving Experimental Reproducibility Through Automated Hybridisation-Based NGS Library Preparation



### Introduction

One of the greatest challenges for laboratories running next generation sequencing (NGS) panels is to ensure high quality, reproducible results whilst sustaining high sample throughput in a timely fashion. Automation enables standardisation of Oxford Gene Technology's (OGT) SureSeq™ NGS library preparation, hybridisation and washing for up to 96 samples in a single batch, improving laboratory productivity. In this application note, an Agilent Bravo® A Automated Liquid Handling Platform was configured to run the SureSeq NGS library preparation protocol. The results demonstrate marked improvement not only in hands-on-time, but also a number of quality metrics, in particular, reduced variance in % on-target reads when compared to manual processing.

### Experimental

Varying amounts of good quality genomic DNA were used as templates. Samples were fragmented using NEBNext® dsDNA Fragmentase® (New England Biolabs Inc. MO348S), libraries were prepared with an Agilent Bravo A Automated Liquid Handling Platform using the SureSeq NGS Library Preparation Kit and hybridised to the SureSeq myPanel™ Custom Myeloid Panel – 49 gene Plus (602017). Resulting libraries were grouped into 16 sample batches and sequenced on an Illumina™ MiSeq (2 x 150 bp v2 cartridge MS-102-2002). Sequencing data was processed using OGT's NGS analysis software Interpret.

### Results

#### DNA Fragmentation

In hybridisation-based target enrichment, genomic DNA is first randomly sheared into 150 – 250 bp fragments by mechanical or enzymatic fragmentation. OGT offers scripts compatible with both mechanical and enzymatic fragmentation workflows, producing highly reproducible results across a range of starting input amounts (Figure 1).

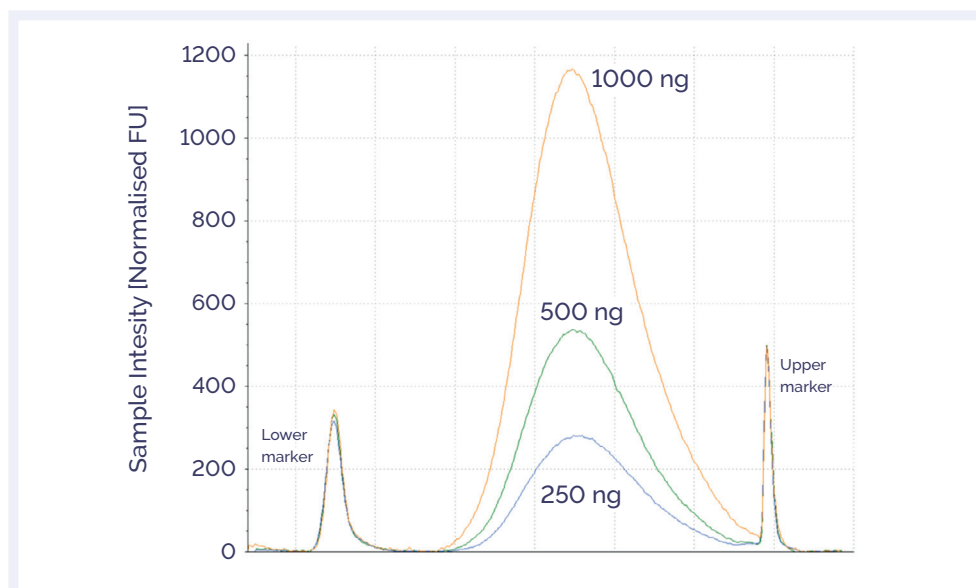


Figure 1: Agilent TapeStation™ trace of enzymatically sheared good quality DNA. Master mix preparation and sample handling, including AMPure® purification, were performed on the Bravo instrument. The fragment sizes generated ranged from 150–250 bp regardless of amount of starting DNA. The recommended fragment size for the SureSeq NGS range is 150–250 bp.

### Reduction in hands-on time

Reducing the labour intensity of any step in any laboratory workflow is highly desirable. As the hybridisation-based enrichment workflow contains routine steps such as reagent/sample transfers and bead-based purifications, it is highly amenable to automation and the benefits that brings compared to manual processing, for example in the reduction in the amount of hands-on time required (Figure 2).

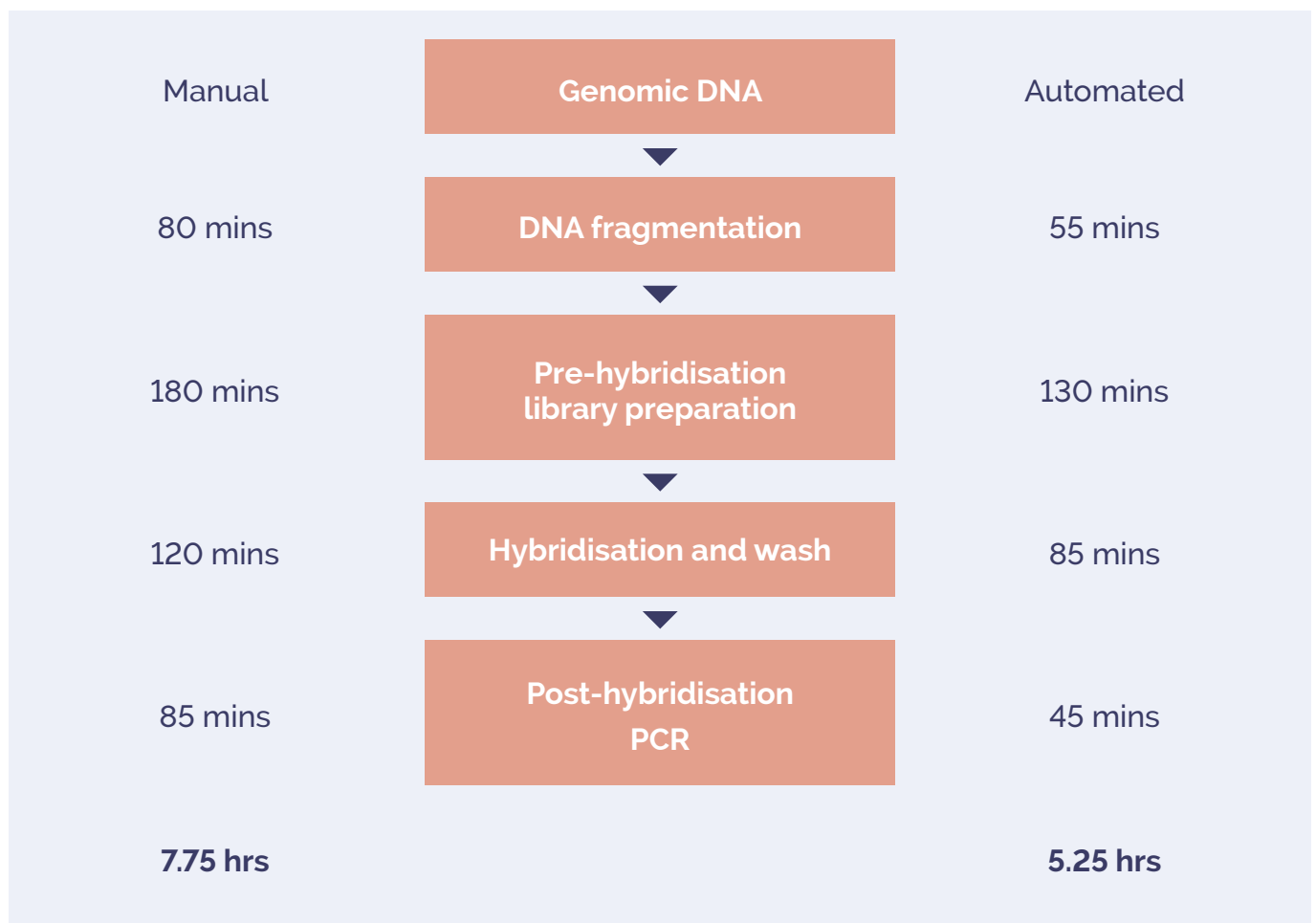


Figure 2: Hands-on time required to process 96 samples using enzymatic fragmentation through to sequence ready libraries.

### Improved reproducibility

NGS hybridisation and washing involves a significant number of pipetting steps which, when performed manually, can introduce high levels of variability. Apart from the requirement for highly skilled personnel, large numbers of manual steps can result in user fatigue, and in some cases errors. Automation offers a solution which can minimise or eliminate these issues, in particular that of increased experimental variability. Samples processed by OGT's scripts are highly reproducible across a range of starting inputs (Figure 3) and show a reduction in variability compared to manually processed samples (Figure 4).

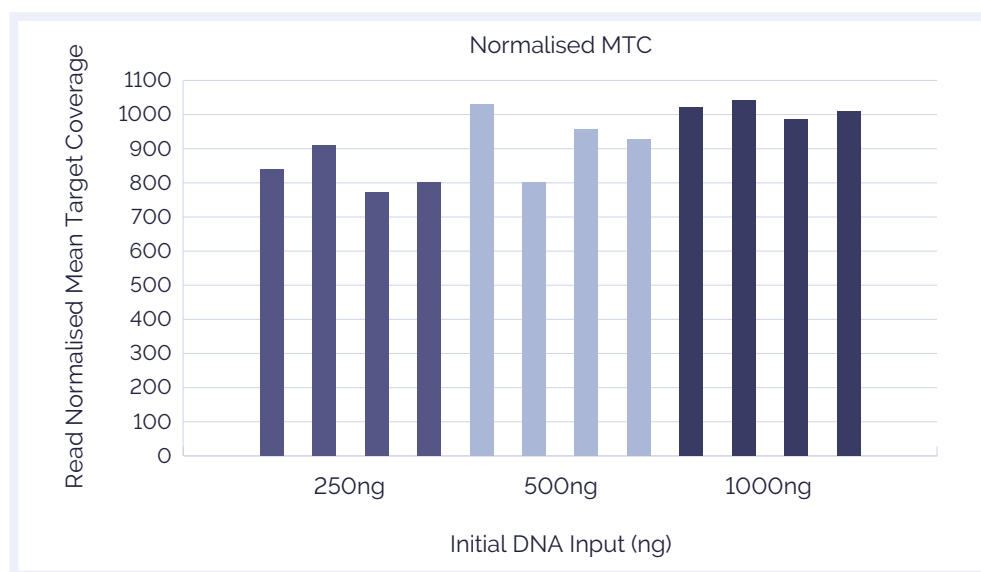


Figure 3: Samples processed on the Agilent Bravo show excellent reproducibility in mean target coverage regardless of the sample input.

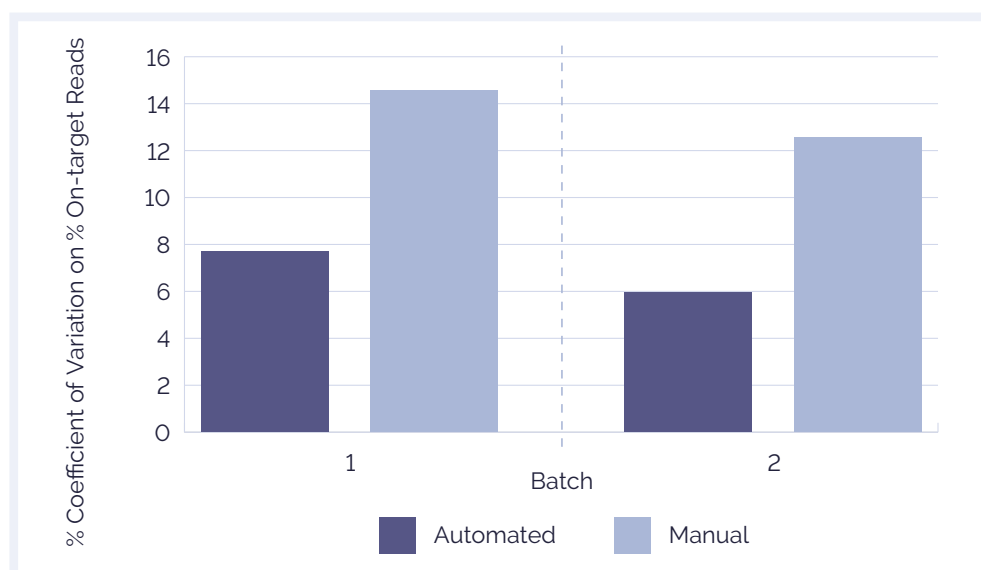


Figure 4: Comparison of a batch of eight samples processed by hand and by automation. The two batches processed by automation show greatly reduced coefficient of variation in % on-target reads compared to two manual batches.

### Compatible with bead-based concentration workflows

Hybridisation-based enrichment protocols often require users to concentrate their unenriched samples using a SpeedVac® or similar vacuum concentration device; this adds time to the total workflow. OGT have developed a bead-based method to concentrate unenriched samples, this alternative workflow is compatible with the Agilent Bravo and delivers highly reproducible data across a range of starting DNA inputs (Figure 5).

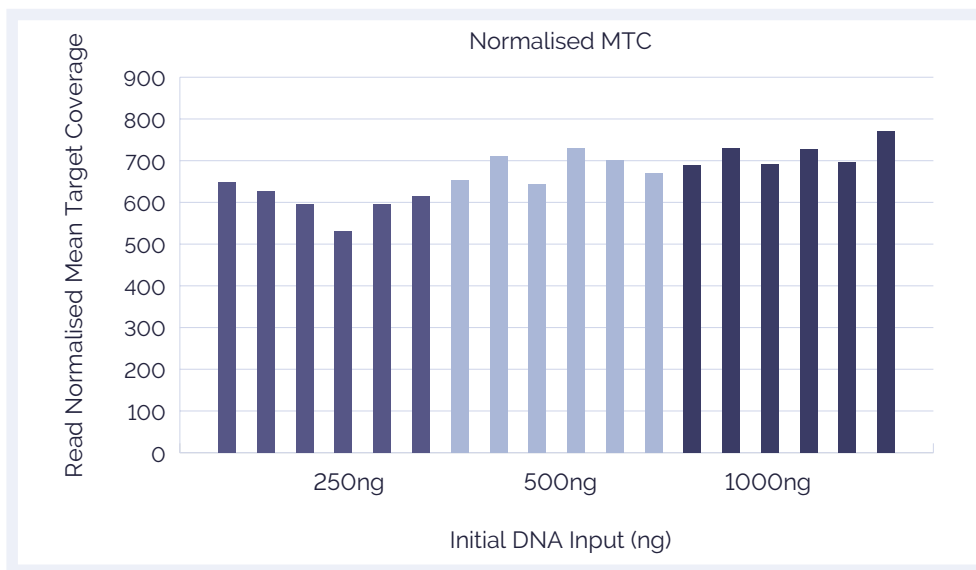
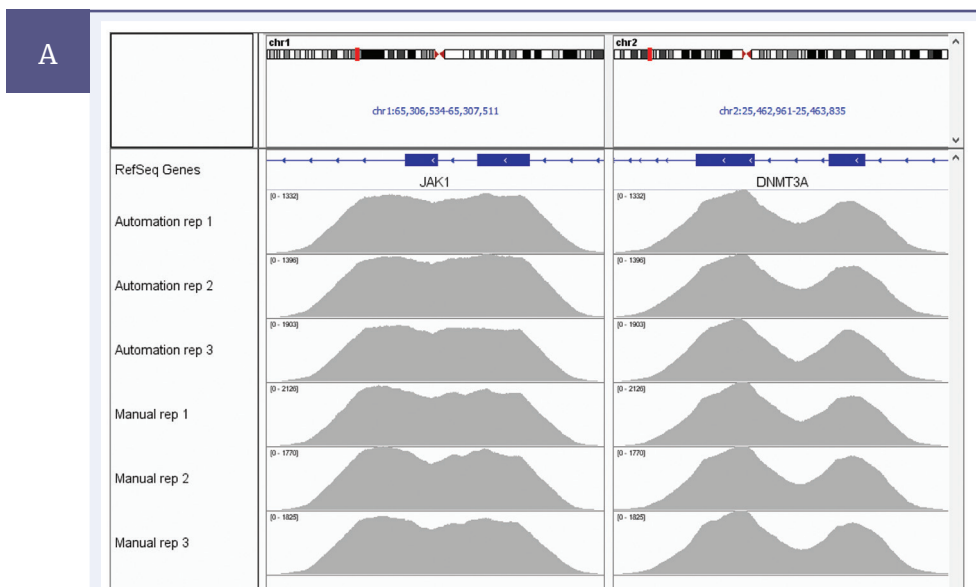


Figure 5: Samples processed using the alternative bead-based concentration workflow.

### High quality sequencing data

Samples processed by the Agilent Bravo maintain and sometimes improve upon the excellent levels of uniformity seen over a variety of targets (Figure 6) including difficult targets such as the *FLT3* ITD region as shown in Figure 7.



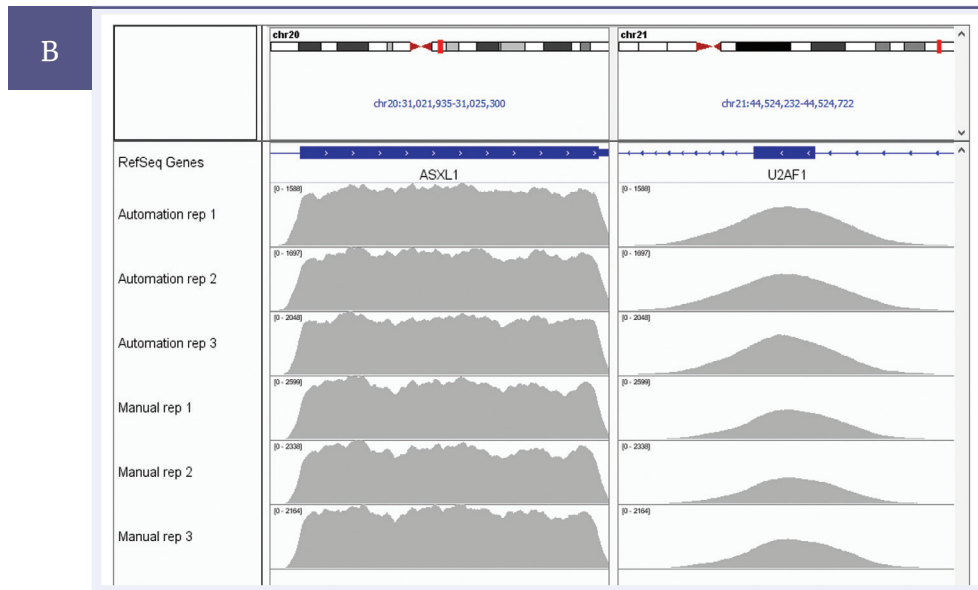


Figure 6: Integrated Genomics Viewer (IGV) image of the coverage over key targets in the SureSeq myPanel Custom Myeloid Panel – 49 gene plus, **A** *JAK1* and *DNMT3A*, **B** *ASXL1* and *U2AF1*.

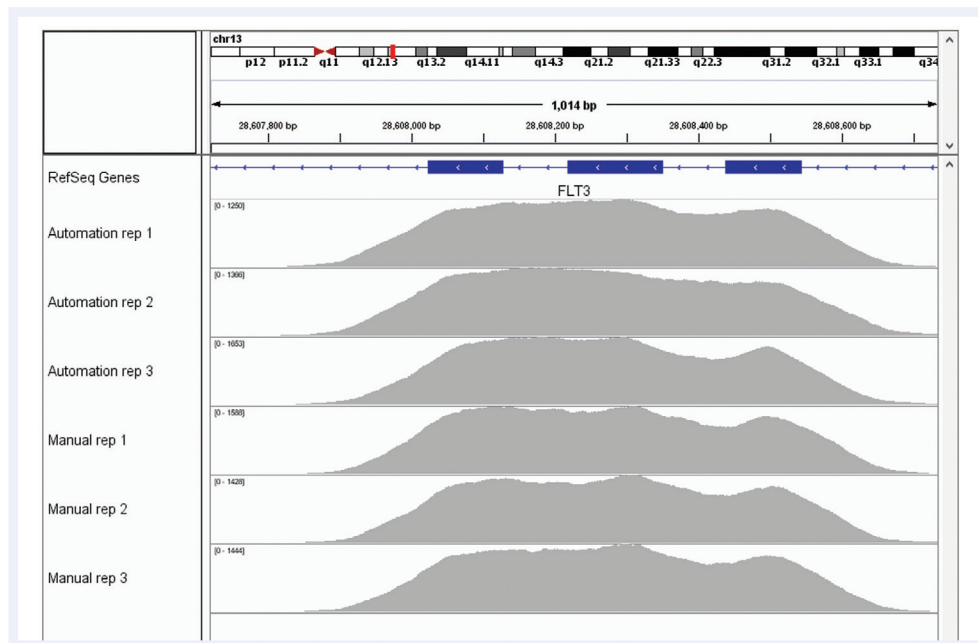


Figure 7: IGV image comparing the uniformity of coverage over *FLT3* exons 13 – 15 in three samples processed by the Agilent Bravo and three manually processed samples, below.

### Conclusions

The Bravo A Automated Liquid Handling Platform was successfully set up to run the SureSeq NGS library preparation workflow. The results showed improved reproducibility with a threefold reduction in the coefficient of variation in % on-target reads of when compared to manual processing, whilst hands on time was reduced by a third.

The ultimate goal of any sequencing assay is to discover all variants present; uniformity of enrichment means that all regions are represented more equally. With this in mind, the excellent reproducibility of mean target coverage and coverage uniformity exhibited by the hybridisation-based enrichment of SureSeq panels in conjunction with the Bravo Platform facilitates the reliable detection of low-frequency somatic variants. With automated preparation allowing the processing of up to 96 samples at any one time, it offers excellent reliability and reproducibility with easy scale up of laboratory throughput.

### Ordering information

**UK +44 (0) 1865 856800**

**US +1 914 467 5285**

**contact@ogt.com**

**ogt.com**



A Sysmex Group Company

**What binds us,  
makes us.**

**Oxford Gene Technology Ltd., Begbroke Science Park, Woodstock Road, Begbroke, Oxfordshire, OX5 1PF, UK**

**SureSeq:** For Research Use Only; Not for Diagnostic Procedures. This document and its contents are © Oxford Gene Technology IP Limited – 2021. All rights reserved. OGT™ and SureSeq™ are trademarks of Oxford Gene Technology IP Limited. The SureSeq NGS Library Preparation Kit was jointly developed between Oxford Gene Technology and Boline Reagents Limited. SpeedVac® is a registered trademark of Thermo Fisher Scientific. Covaris™ is a registered trademark of Covaris, Inc., Bravo® and TapeStation™ are registered trademarks of Agilent Technologies, Inc., NEBNext® and Fragmentase® (New England BioLabs, Inc.) are registered trademarks of New England Biolabs, Inc. AMPure® is a registered trademark of Beckman Coulter, Inc.