

## myProbes Passport (RUO Probes)

<b>MPP2900</b>	<b>NTRK3/ETV6 Dual Fusion Probe</b> <i>For Research Use Only</i>
----------------	---



**Nucleotide Locations** (If requesting more than one probe please specify individual probe colours)

	<b>Chromosome Locus</b>	<b>Start</b>	<b>End</b>	<b>Colour</b>
1	15q25.3	88171760	88394033	Red
2	15q25.3	88328543	88494404	Red
3	15q25.3	88621692	88794005	Red
4	15q25.3	88745702	88907307	Red
5	15q25.3	88864972	89072568	Red
6	12p13.2	11526867	11694579	Green
7	12p13.2	11692825	11801621	Green
8	12p13.2	11818561	11946246	Green
9	12p13.2	11924600	12041067	Green
10	12p13.2	12055716	12223558	Green
11	12p13.2	12228893	12420408	Green

Database Used and Version

GRCh37 (hg19)

## Probe Map

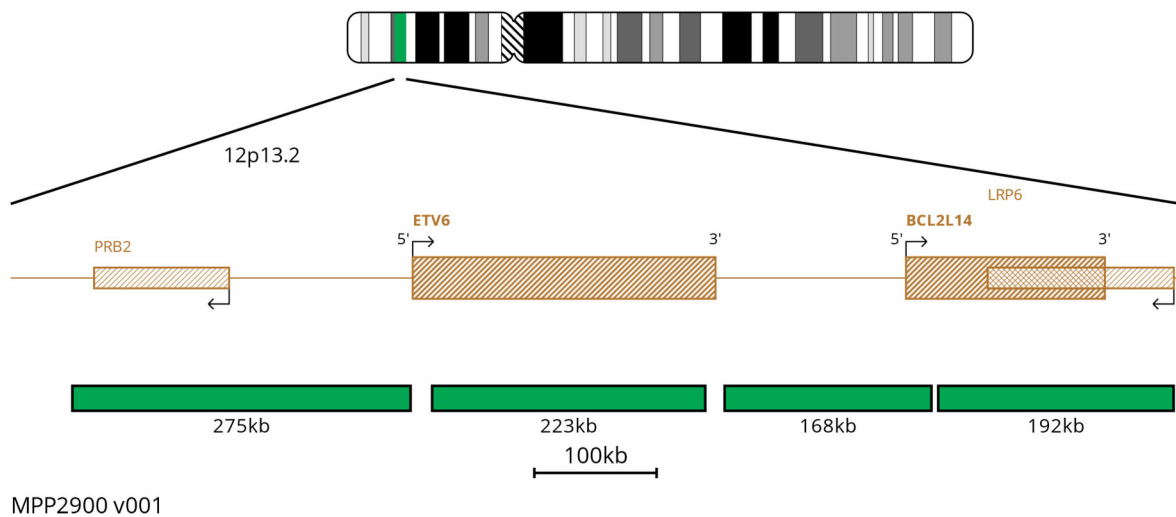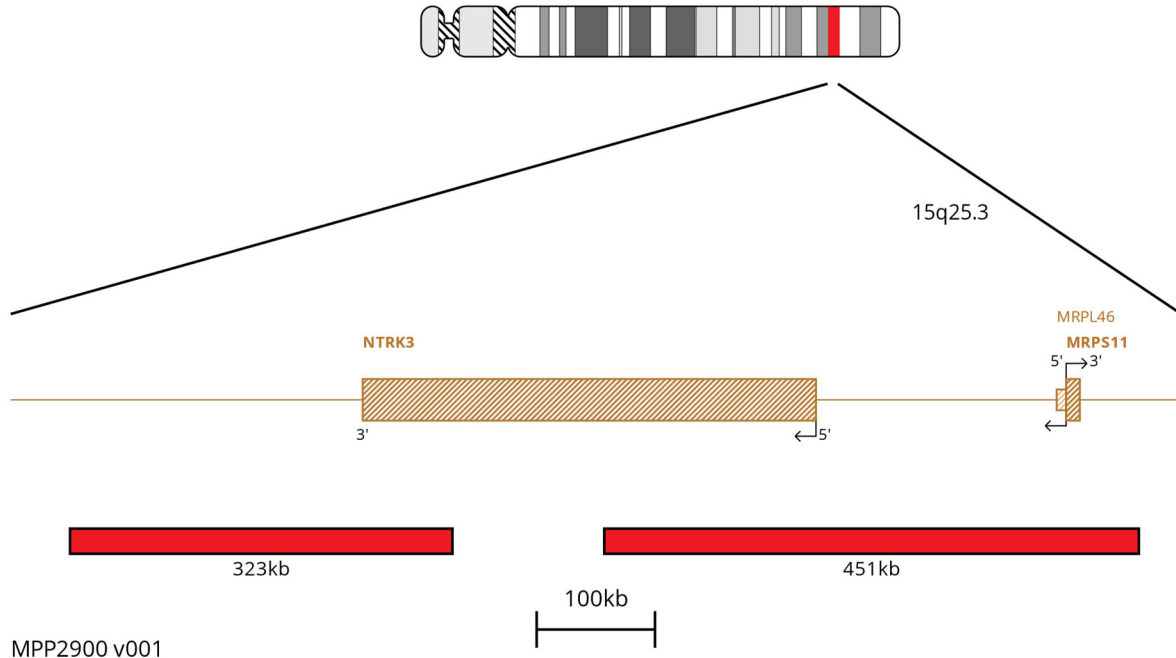

## Materials Provided

- Probe:** 100µl per vial  
 The probe is provided pre-mixed in hybridisation solution (Formamide; Dextran Sulphate; SSC) and is ready to use. It is directly labelled with a red fluorophore, which lies in the Texas Red Spectrum and a green fluorophore, which lies in the FITC Spectrum. The probe has been batch released after QC testing on Formalin Fixed Paraffin Embedded Sample.
- Counterstain:** 150µl per vial. The counterstain is DAPI antifade (ES: 0.125µg/ml DAPI (4,6-diamidino-2-phenylindole)).

## Warnings and Precautions

1. For research use only. Not for use in diagnostic procedures.
2. For laboratory professional use only.
3. Probe mixtures contain formamide, which is a teratogen; do not breathe fumes or allow skin contact. Handle with care; wear gloves and a lab coat.
4. DAPI is a potential carcinogen. Handle with care; wear gloves and a lab coat.
5. Follow local disposal regulations for your location along with recommendations in the Safety Data Sheet to determine the safe disposal of this product. This also applies to damaged test kit contents.
6. Dispose of all used reagents and any other contaminated disposable materials following procedures for infectious or potentially infectious waste. It is the responsibility of each laboratory to handle solid and liquid waste according to their nature and degree of hazardousness and to treat and dispose of them (or have them treated and disposed of) in accordance with any applicable regulations.
7. Operators must be capable of distinguishing the colours red, blue, and green.
8. The probe should not be diluted or mixed with other probes.
9. All products should be validated before use.
10. Internal controls should be carried out by using unaffected cell populations in testing samples.
11. Custom probes are specifically developed for individual customers' RESEARCH USE ONLY (RUO) requirements and not with the intention of being used for in vitro diagnostic examination. Therefore, prior to any use of these probes, users should review the design of such probes to confirm they are suitable for their requirements.

## Storage and Handling

1. Store the probe between -25°C to -15°C.
2. Based on the stability established for other substantially equivalent CytoCell probes, this myProbes product should be stable for 2 years post manufacture date when stored as indicated on the label.
3. Store the probe and counterstain vials in the dark. Ensure that exposure of the probe and counterstain to laboratory lights is limited at all times.

## Known Cross-Reactivity

No cross hybridisation observed at QC release on the sample type specified.

## Customer Support

Please contact the CytoCell Technical Support Department or email [probes@cytoCELL.com](mailto:probes@cytoCELL.com).

## Patents and Trademarks

myProbes and CytoCell are registered trademarks of CytoCell Ltd.

For Research Use Only. Not for use in diagnostic procedures.

This product contains technology licensed from Life Technologies Corporation that is available for human diagnostics or life science research use only.



A Sysmex Group Company

**CytoCell Ltd.**

Oxford Gene Technology  
418 Cambridge Science Park,  
Milton Road,  
Cambridge, CB4 0PZ, UK  
T: +44(0)1223 294048  
F: +44(0)1223 294986  
E: [probes@cytoCell.com](mailto:probes@cytoCell.com)  
W: [www.ogt.com](http://www.ogt.com)

